




operating table SU-02

Purpose

Designed for procedures in the range of general surgery, vascular surgery, cardiosurgery, neurosurgery, urology, proctology, gynaecology, laryngology, ophthalmology, endoscopy, laparoscopy, traumatology, oncology, dentistry, plastic surgery and others.

Full range of specialist equipment is introduced inside "Equipment for operating tables" catalogue.

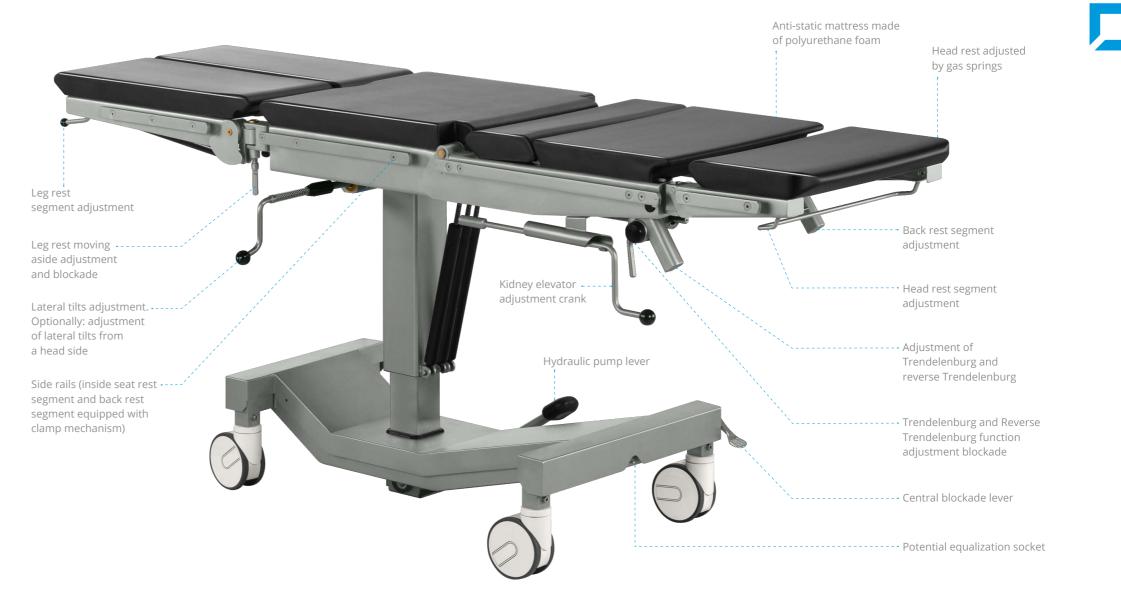

Standard structure

- The table is constructed of acid proof, stainless steel with mat table surface
- The mobile table base is operated by a foot lever (In the version equipped with base made in shape of T letter locking of the table is done with the use of four extended foots. In the operating table version equipped with base made in shape of H letter here is a central locking system equipped with directional wheel)
- Operating table surfaces are easy to clean and resistant to disinfecting agents
- There is a possibility to change head rest segment with leg rest segment
- Possible use of plate leg rest segment type

Functions and drives

- Table top height adjustment with use of hydraulic foot pump
- Back rest segment angle inclination, leg rest segment (regarding four and five segments table top), head rest segment, Trendelenburg and Reverse Trendelenburg functions are supported by gas springs with blockade
- Kidney elevator (in case of five and six segment table top adjusted manually with use of crank placed at patient's head side)
- Lateral tilts are executed with the use of a crank, located on the side of the table under the seat segment or optionally by a crank located at the patient head side
- In case of a six segments table top leg rest segments are manually adjusted thanks to gear transmission
- Leg rest segments are manually aside adjusted

SU-02 operating table bases

- Operating table SU-02 equipped with base made of T letter shape
- Operating table SU-02 equipped with base made of H letter shape
- Operating table SU-02 without of lateral tilts equipped with base made of H letter shape

Operating table SU-02 equipped with base made of H letter shape



Operating table SU-02 equipped with base made of T letter shape



Operating table SU-02 without of lateral tilts equipped with base made of H letter shape



Standard functions and solutions

Co-operation with C arm

Possibility of exchanging leg rest segments with head rest segment



X-rays permeability



Range of X-ray permeability with use of the RTG tray



Possibility mattresses removal

Table tops SU-02:

Four segments table top with non-divided leg rest:

total length: 2035 mm

- total width: 550 mm
- leg segment length: 615 mm

Four segments table top:


- total length: 2035 mm
- 550 mm • total width:
- leg segment length: 615 mm

Five segments table top:

- total length: 2035 mm
- total width: 550 mm
- leg segment length: 615 mm

Six segments table top:

- 2160 mm total length:
- 550 mm • total width:
- leg segment length: 750 mm

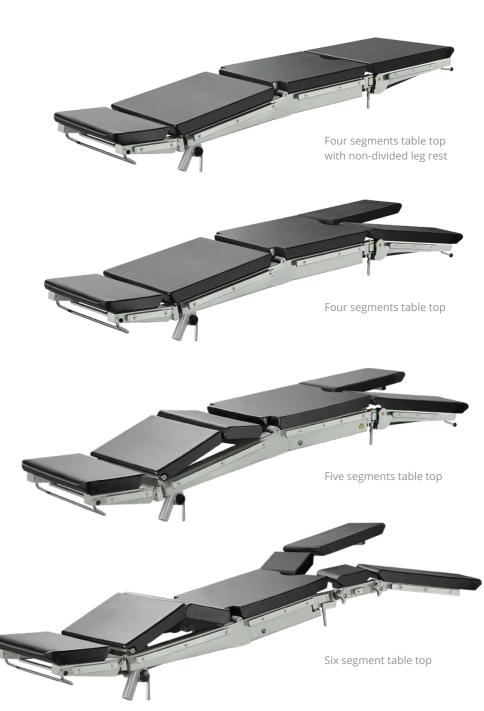

Trendelenburg and reverse Tredelenburg position







Lateral tilts



SU-02 execution options

- Operating table SU-02 equipped with base made of T letter shape
- Operating table SU-02 equipped with base made of H letter shape
- Operating table SU-02 without of lateral tilts equipped with base made of H letter shape
- Operating table SU-02 equipped with lateral tilts located from patient head side
- Operating table SU-02 equipped with plate leg rest segment



equipped with lateral tilts located from patient head side



Operating table SU-02 equipped with plate leg rest segment



Operating table SU-02 without of lateral tilts equipped with base made of H letter shape

Examples of use



Patient positioning for laparoscopy position (moved aside and lowering down leg rest segments allow for anatomic patient positioning for treatment)



Gynaecology position equipped with gynaecology attachment WS-28.7 and knee crutches WS-05.5

Orthopaedic attachment

Orthopaedic attachment made of stainless acid proof steel is fixed for table top. There is possible to carry out together with table top longitudinal and lateral shifts and height adjustment. It's possible patient monitoring by C arm and ensures for surgeon comfortable access for a patient









Kidney elevator (in case of five and six segment table top adjusted manually with use of crank placed at patient's head side)



Seating position with use of specialist head rest WS-21.5

Technica data	with base made of T letter shape	with base made of H letter shape	with base made of H letter shape (without of lateral tilts)	FAI Medical Sc
Total table length equipped with 4 or 5 segment table top	2035 mm	2035 mm	2035 mm	Famed Żywiec ul. Fa 34-300 Żywie Hea tel.: +48 33 fax: +48 33 Export Dep tel.: +48 33 tel.: +48 33 tel.: +48 33 fax: +48 33
Total table length equipped with 6-segment table top	2160 mm	2160 mm	2160 mm	
Total table top width	550 mm	550 mm	550 mm	
Mattress width x thickness	500 x 50 mm	500 x 50 mm	500 x 50 mm	
Table top height adjustment	725 to 1025 mm	745 to 1045 mm	725 to 1025 mm	
Trendelenburg position	30°	30°	30°	
Reverse Trendelenburg position	30°	30°	30°	
Lateral tilts	±25°	±25°	Lack of	
Back rest segment adjustment	-45° to +85°	-45° to +85°	-45° to +85°	
Head rest adjustment	-50° to +55°	-50° to +55°	-50° to +55°	expo.famed@fam www.fam
Kidney elevator height adjustment (5 and 6 segments table top)	150 mm	150 mm	150 mm	
Leg rests Inclination angle adjustment (4 and 5 table top segment)	-90° to +25°	-90° to +25°	-90° to +25°	
Leg rests Inclination angle adjustment (6 table top segments)	-90° to +90°	-90° to +90°	-90° to +90°	
Leg rest segment Inclination aside angle adjustment	180°	180°	180°	
Permissible load (leg rest segment fixed in seat segment / leg rest segment in back rest segment	200 / 135 kg	200 / 135 kg	200 / 135 kg	
Table weight	160 kg	105 kg	105 kg	

Operating table functions performed by foot pump

÷

725 to 1025 mm (SU-02 with base made of T and H without of lateral tilts) 745 to 1045 mm (SU-02 with base made of H letter shape)



CE


Operating table functions performed manually

+25°



5



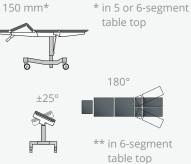

ec Sp. z o.o. abryczna 1 /iec, Poland

ead Office: 3 866 62 00 3 475 58 90

partment:

3 866 62 58 3 866 62 65 3 866 62 66 3 861 46 78

med.com.pl med.com.pl

The product meets the requirements of the European Directive MDD 93/42/EEC concerning safety of medical devices



Famed Żywiec reserves the right to make modifications in the products in the course of technical progress. All illustrations and photos, used in this material, are for illustration and may not accurately depict the actual product.